

## CORN DETECTION FROM FOOD

### REAL TIME PCR KIT

Cat. No: 313R-10-01

#### PRODUCT DESCRIPTION

According to the European Commission Directive 2002/86/EC, food ingredients have to be declared. SNP Animal Species Detection for Corn Real Time PCR Kit is used for the detection of Corn (*Zea mays*) DNA in raw and processed food. The Kit provides high sensitivity with a limit of detection (LOD) below the officially required level of 0.1%.

#### PRINCIPLE OF THE SYSTEM

During the PCR reaction, the DNA polymerase cleaves the probe at the 5' end and separates the reporter dye from the quencher dye only when the probe hybridizes perfectly to the target DNA. This cleavage results in the fluorescent signal which is monitored by Real-Time PCR detection system. An increase in the fluorescent signal (CT) is proportional to the amount of the specific PCR product.

#### PRODUCT SPECIFICATION

Each isolated DNA should be tested with Corn Real Time master mix. The system provides reagents in a "Ready-to-Use" master mix format which has been specifically adapted. The test system is designed for use with sequence specific primers and probe for mitochondrial DNA of Corn.

The fluorescence of species analysis is FAM. Also mastermix contains an internal control labelled with HEX/JOE dye for check PCR reaction.

#### SYSTEM CONTENTS

Reagents	20 rxns
• Corn Real Time PCR Master Mix	400 µl
• Corn Positive Control DNA	30 µl
• Corn Negative Control DNA*	30 µl

(\*): Soybean DNA for check Cross-Reaction.

#### STORAGE

- All reagents should be stored at – 20 °C and dark.
- All reagents can be used until the expiration date on the box label.
- Repeated thawing and freezing (>3X) should be avoided, as this may reduce the sensitivity of the assay.

#### DNA EXTRACTION

Corn DNA can be isolated by using any DNA isolation protocol from plant. Our system optimized according to CTAB DNA isolation protocol from plant.

#### PROCEDURE

- Before starting work, mix the mastermixes gently by pipetting
- For each sample, pipet **20 µl mastermix\*** with micropipets of sterile filter tips to each optical white strips or tubes.
- Add **5 µl DNA** into each tube.
- Run with the programme shown below.

\*Master mixes include *HotStart Taq DNA Polymerase*.

#### PCR PROGRAMME

95 °C	3 Min.	Holding
95 °C	15 Sec.	30 Cycles
60 °C	1 Min.	

Fluorescent dyes are FAM and HEX/JOE

#### If you use:

- ABI Prism® system, please choose **"none"** as passive reference and quencher.
- Mic qPCR Cycler, please adjust gain settings, **"Green Auto Gain"** to **20** and **"Yellow Auto Gain"** to **10**.


#### This system can use with:

Bio-Rad CFX96  
ABI Prism® 7500/7500 Fast  
Roche LightCycler® 480 System  
Rotor Gene Q  
Mic qPCR Cycler



## DATA ANALYSIS

After the run is completed data are analysed using the software with HEX (JOE) and FAM dyes. The below results were studied with ABI 7500.

Internal control amplification plots must be seen in all wells and has been labelled with HEX/JOE dye. The CT value of internal controls shows  $22 \leq X \leq 27$ .

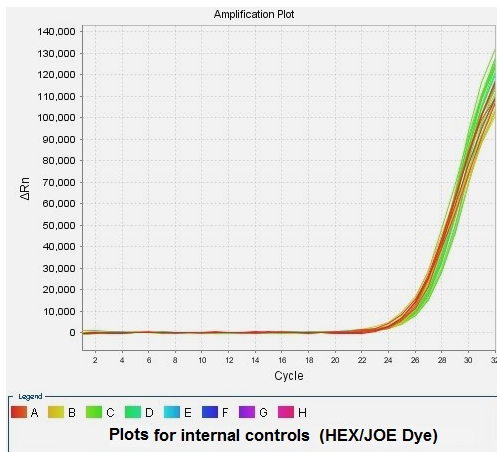

Corn DNA can be analyzed with FAM dye. The CT value of positive sample shows  $17 \leq X \leq 29$ .

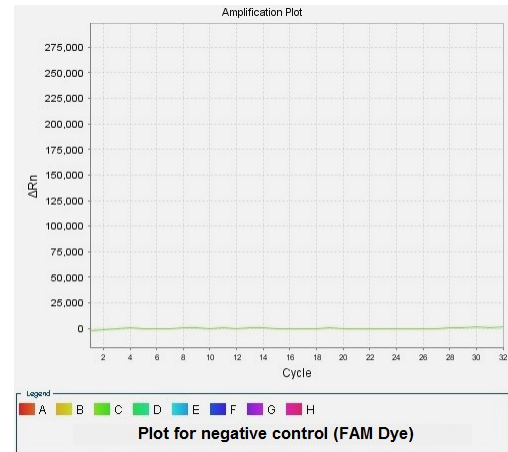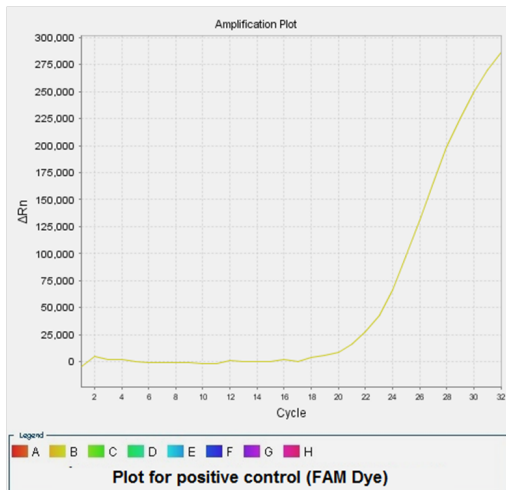

## TROUBLE SHOOTING

### If internal control doesn't work,

- Absence of DNA
- Sample is containing DNA inhibitor(s)

### If plots start late,

Compare positive control and sample. If there is no problem in positive control,

- DNA quality is not good.
- The amount of DNA is not enough.

Please contact us for your questions. [tech@snp.com.tr](mailto:tech@snp.com.tr)

## CAUTIONS

- All reagents should be stored at suitable conditions.
- Do not use the PCR mastermixes forgotten at room temperature.
- Thaw PCR mastermix at room temperature and slowly mix by inverting before use.
- Shelf-life of PCR mastermix is 12 months. Please check the manufacturing data before use.
- Only use in vitro diagnostics.